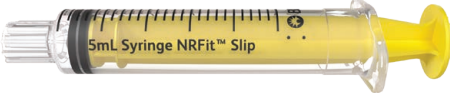


BD® Syringes NRFit™ Slip and NRFit™ Lok

Tip sheet

What are the changes?

Changes	BD Plastipak™ Syringes	BD® Syringes NRFit™ Lok and NRFit™ Slip	
Connector design	Universal Luer connector that can connect to any Luer general use product.	• Neuraxial-specific connector design, incompatible with other applications. • Reduce drug waste compared to BD Plastipak™ Syringes.	
Lock syringes		Require an extra half-turn rotation on average for a secure connection compared to Luer lock devices, for added confidence in the security of the connection.	
Slip syringes	No collar. BD Plastipak™ Luer Slip Tip Syringe 	• Additional collar to slip syringes. • Can still be connected to a NRFit™ needle with a push connection.	BD® Syringe NRFit™ Slip 
Syringe labeling distinguishing slip and lock syringes	No labeling on the product itself denoting slip or lock, but the syringe tips are easily distinguishable. Lok:  Slip: 	New “Lok” visual indicator to help distinguish BD NRFit™ Lok Vs BD NRFit™ Slip. BD Syringe NRFit™Lok:  	BD Syringe NRFit™ Slip: 
Syringe colour coding	Not coloured.	NRFit™ packaging label and yellow component.	BD® Syringe NRFit™ Slip with yellow plunger rod 

1030 Eskdale Road, Winnersh Triangle, Wokingham, RG41 5TS

[bd.com/en-uk](https://www.bd.com/en-uk)

